
M.Sc.Dr. Khagani Huseynov

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He graduated from Azerbaijan Medical University Faculty of Dentistry in 2001. In 2006, he completed the expertise course on Maxillofacial Surgery and Implantology in Post-Diploma Training Academy, St. Petersburg University. In 2011, he received his master degree in laser dentistry at Achen University, Germany. In 2012, he completed the surgical courses of Prof. Dr. Ateş Parlar. In 2015 he graduated from the education program of Implantology in Frankfurt am Main University, Germany. He is also the head of the Azerbaijan Esthetics Dentistry Academic Association). He has been various organising congresses and giving seminars and courses in his country since 2008.

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
*Patient Information:
Patient is a male with poor mouth hygiene and was born on August 3, 1969.
Patient History:
The patient consulted to the clinic with pain and loose teeth.*

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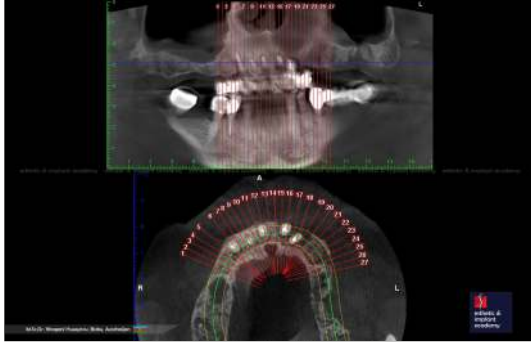
Pre-operation CBCT(Cone Beam Computed Tomography) results have shown maxillary teeth are inflammatory. The teeth have fourth degree loose and decided for extraction.

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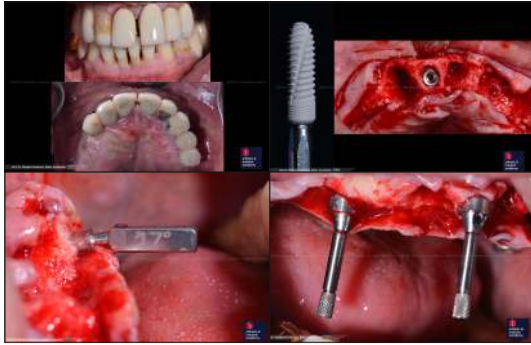The old metal porcelain teeth in mouth before surgery. After the extraction, flep has been removed, socket has been curetted from infected tissues. In operation, angular pins used to provide parallelism of multi unit abutments. Nucleoss T6 4.1/14 mm implant has been used by mounting 1 mm depth.

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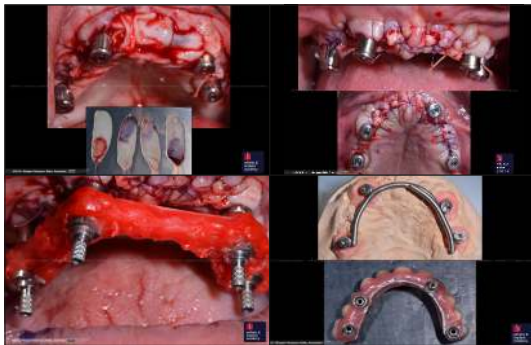L-PRF has used to accelerate recovery period in extraction sockets. After inserting healing caps, 3/0 soluble P.G.A. suture has used to close the scar. The pictures of the temporary prosthetics.

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After four months, the size is measured to prepare the permanent teeth. The pictures of the process of preparing the base platform, plastic caps before sending metal casting.

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Esthetic result and satisfied patient.

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X-ray scan of final situation after four months. Implants are healthy, no bone resorption at all, on the quite contrary there is bone incresement due to platform switching. Prosthesis adaptation is perfect!